



HFMED®

**YD01-4** | Emergency Cold light operating lamp

**Instructions for use:**

For basic health infrastructures; in support of surgical interventions in operating theatres, delivery rooms, dressing rooms and emergency departments.

After each intervention the handle should be disinfected, cleaned and sterilized by steam sterilizer.

Never touch a halogen bulb with bare hands.

Must be used and maintained properly, according manufacturer's instructions.



**WHY CHOOSE US**



**OUR ADVANTAGE**

The company keeps to the concept of being oriented by science and technology and taking quality as life. Now the company has passed the certification of ISO9001:2008, ISO13485:2003, CE and SGS.



**MANAGEMENT**

"Quality first, Credit first, Service first" is our tent. Huifeng are looking forward to sincere cooperation with you. Let's advance and pursue a wonderful future by mutual efforts.



**SALES MARKET**

all has its domestic sales amount coming out at the top and some products have ready sale in Southeast Asia. Its products can be found in all



Official QR code



HFMED®

YOUR IDEAL SUPGICAL LIGHT SUPPLIER  
**SHANGHAI HUIFENG MEDICAL INSTRUMENT CO., LTD.**

Address: NO.885, Qiuxing Road, Nicheng Town, Pudong New Area, Shanghai, China.  
 Tel: +86-021-58084992  
 Mob: +86-13918165133  
 Email: lisa@huifeng-co.com  
 web: http://en.huifeng-co.com


HFMED®

CE




**YD01-4**

Emergency Cold light operating lamp

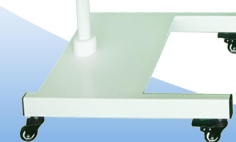

**HFMED**

SHANGHAI HUIFENG MEDICAL INSTRUMENT CO.,LTD

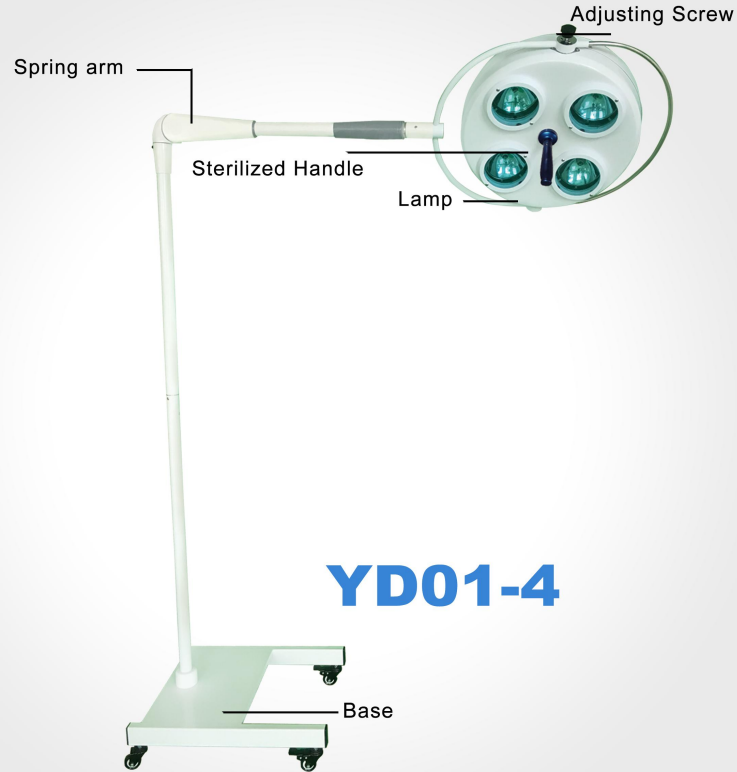

## YD01-4

### Product Description

Device used for providing a specialized source of light for illumination of a site of medical intervention. It provides a high intensity, high colour rendering field of light that minimizes shadows and the emission of heat. Adjustable Operating field diameter with focus



Automatic switch from mains to batteries in case of power failure  
With battery status indicator  
Removable autoclavable handle  
Power supply: The supplier has to supply the right voltage, plug and socket type that fit the country it is delivering the goods to. The country of destination would be stated in the UNFPA's



Mobile on heavy sturdy stand,  
On 4 antistatic ball-bearing swivel castors, 2 with breaks  
Spring balanced articulating arm  
Low centre of gravity for optimal stability and reach



### Main Technical Data

YD01-4	Shadowless Operating lamp
Illuminance	≥50,000Lux
Colour temperature	3700±500K
Import power	150VA
Illuminance depth	≥600mm
Total irradiance	547W/ m <sup>2</sup>
Size of light field	160-220mm
Lamp head diameter	500mm
Rated power of bulb	25W/24v
Power supply voltage	220V/50Hz
Bulb	cold light
Accessories/ spare parts/consumables	
	1 x Set of spare halogen bulbs
	1 x Spare handle
	1 x Set of spare fuses